



## SET-UP INSTRUCTIONS

- ▶ Connect the **VT UNIT** 18-inch air hose to the air (AIR) filter trap, and the **VT UNIT** 18-inch oxygen hose to the oxygen (O<sub>2</sub>) filter trap on the Precision Flow.
- ▶ Attach second 18-inch air/ O<sub>2</sub> hose to the **TANK** NIST connector on the air/O<sub>2</sub> manifold.
- ▶ Attach the 10-foot air/ O<sub>2</sub> hose to the **WALL** NIST connector on the air/ O<sub>2</sub> manifold.

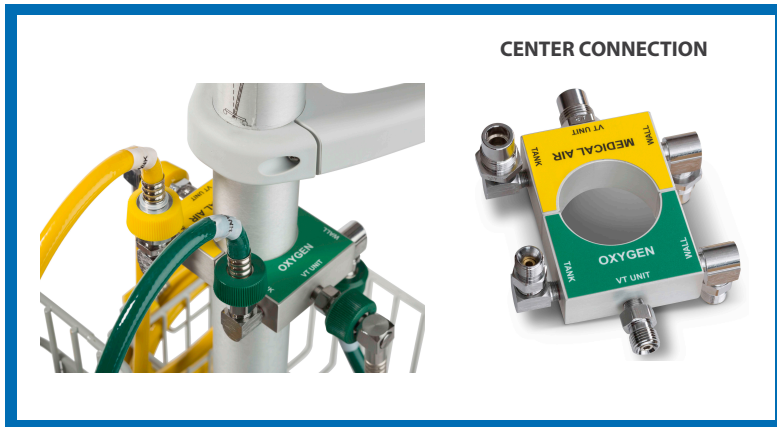

# VAPOTHERM Transfer Unit

## QUICK REFERENCE GUIDE

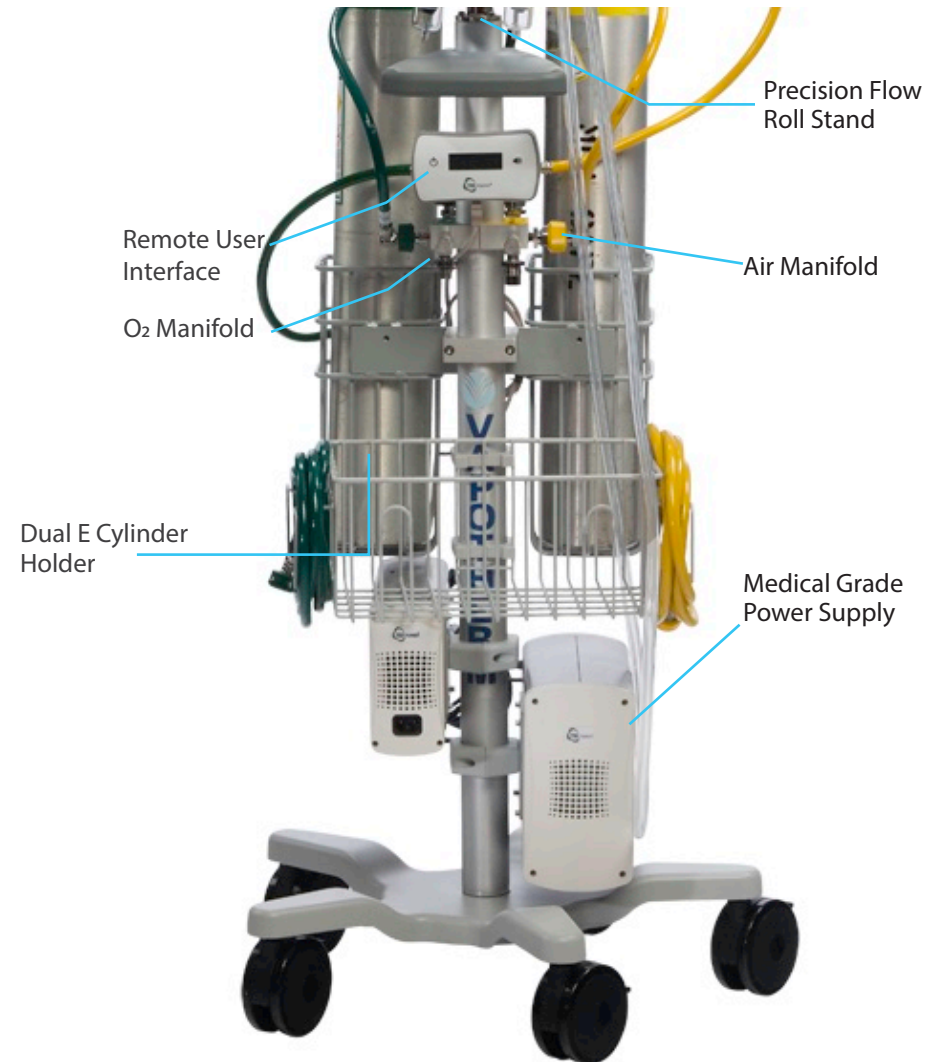

## SWITCHING GAS AND POWER SOURCES

1. Confirm the Precision Flow therapy settings.
2. Open the oxygen and air E cylinders and confirm that the tanks contain adequate gas supply.  
**Warning: Do not attempt to transfer a patient with  $\leq 400$ psi in either tank.**
3. Disconnect the oxygen and air hoses from the wall gas supply.
4. Unplug the power cord from the Astrodyne-TDI power supply and store securely on the unit.
5. Transfer patient to desired location within the hospital.
6. At destination:
  - Connect the oxygen and air hoses to the wall gas supply.
  - Close the oxygen and air E cylinders.
  - Plug the Astrodyne-TDI power supply into a hospital-grade outlet.



**VAPOTHERM, INC.**  
100 Domain Drive  
Exeter, NH 03833  
T: 603-658-0011  
USA

**Technical Support**  
T: +1 (603) 658-5121  
International  
855 557 8276 Domestic  
ts@vtherm.com

May be patented. [www.vapotherm.com/patents](http://www.vapotherm.com/patents)

3101254 Rev. A

CE 0297

EC REP

RMS-UK Ltd.  
28 Trinity Road  
Nailsea, North Somerset BS48 4NU  
United Kingdom



This guide provides you with basic instructions on how to set up and operate the Vapotherm Transfer Unit. The Vapotherm Transfer Unit should only be used with the Vapotherm Precision Flow.

## VTU MOBILE THERAPY DELIVERY TIMES

### LOW FLOW DISPOSABLE PATIENT CIRCUIT

(Duration of use blending from E-size oxygen and E-size air cylinders; times shown in minutes).

TOTAL FLOW L/ MIN	% OXYGEN								
	21%	30%	35%	50%	60%	70%	80%	90%	100%
1	560	632	681	885	1106	903	750	641	560
2	280	316	340	442	553	451	375	321	280
3	187	211	226	295	369	301	250	214	187
4	140	158	170	221	277	226	187	160	140
5	112	126	136	177	221	181	150	128	112
6	93	105	113	147	184	150	125	107	93
7	80	90	97	126	158	129	107	92	80
8	70	79	85	111	138	113	94	80	70

### CLINICAL USE

When transferring a patient between care areas you may swap the disposable patient circuit from the VTU to a stationary Precision Flow unit.

### CARE AREAS

The VTU may be used to transfer patients within acute-care medical environments for which the Precision Flow is cleared.

**Warning: These mobile run time estimates are based on use of a 2000-psi E cylinder. Actual performance may vary depending on the amount of gas in the cylinders.**

### HIGH FLOW DISPOSABLE PATIENT CIRCUIT

(Duration of use blending from E-size oxygen and E-size air cylinders; times shown in minutes)

TOTAL FLOW L/MIN	% OXYGEN								
	21%	30%	35%	50%	60%	70%	80%	90%	100%
5	112	126	136	177	221	181	150	128	112
6	93	105	113	147	184	150	125	107	93
7	80	90	97	126	158	129	107	92	80
8	70	79	85	111	138	113	94	80	70
9	62	70	76	98	123	100	83	71	62
10	56	63	68	88	111	90	75	64	56
15	37	42	45	59	74	60	50	43	37
20	28	32	34	44	55	45	37	32	28
25	22	25	27	35	44	36	30	26	22
30	19	21	23	29	37	30	25	21	19
40	14	16	17	22	28	23	19	16	14

The Vapotherm Transfer Unit (VTU) allows the mobile delivery of high-flow humidified respiratory therapy within a hospital environment. With a fully charged battery module and filled gas supply tanks, the duration of available mobile therapy depends on 2 factors:

- **Tank gas supply:** mobile run time will vary depending on the O<sub>2</sub>/air blend and flow rate settings on the Precision Flow.
- **Battery life:** Up to 1 hour from a full charge. Do not leave the battery charging for a prolonged period of time to avoid adverse battery performance by overcharging. Leaving the batteries fully discharged for long periods of time DOES adversely affect battery performance. For long term storage, remove the cables from the battery after it has been fully charged. The battery should be charged at least once every 4 months.